

### MODEL GL115 PATHOLOGY WORKSTATION



### MODEL GL115

The Model GL 115 has the same useful features as the Model GL 110 with the added bonus of a larger sink basin. The larger sink basin allows the unit to be equipped with an optional disposal unit.

#### Standard Design Features:

- Available in two sizes: 48 inches long and 60 inches long
- All stainless steel construction
- Integral sink with mixing faucet
- Top mounted fluorescent light fixture
- Two-incandescent lights with 3x magnifier mounted on flexible arms
- Formalin container with spigot
- C-fold paper towel holder, stainless steel
- Magnetic instrument holder
- Polyethylene dissecting board



**MODEL GL115 DRAWING, OPTIONAL FEATURES & FILTERS**


**MODEL GL115**


**Optional Accessories:**

- Heavy-Duty commercial disposal (1/2 hp, 3/4 hp, 1 hp, 2 hp)
- Exhaust duct connection to building ventilation system
- Digital Platform Scale, 2 kg capacity, economy or stainless grades available
- Dictation Equipment with foot control or voice operation
- Camera mount facilities
- Dissecting Area Rinse Assembly
- Eyewash Assembly
- Model GL 120 or Model GL 125 Pathology Station Bases
- Flammable Storage Cabinet

**Replacement Filters:**

- Model GL 112 Activated Charcoal Filter for odorous fumes
- Model GL 113 Potassium Permanganate for formaldehyde fumes

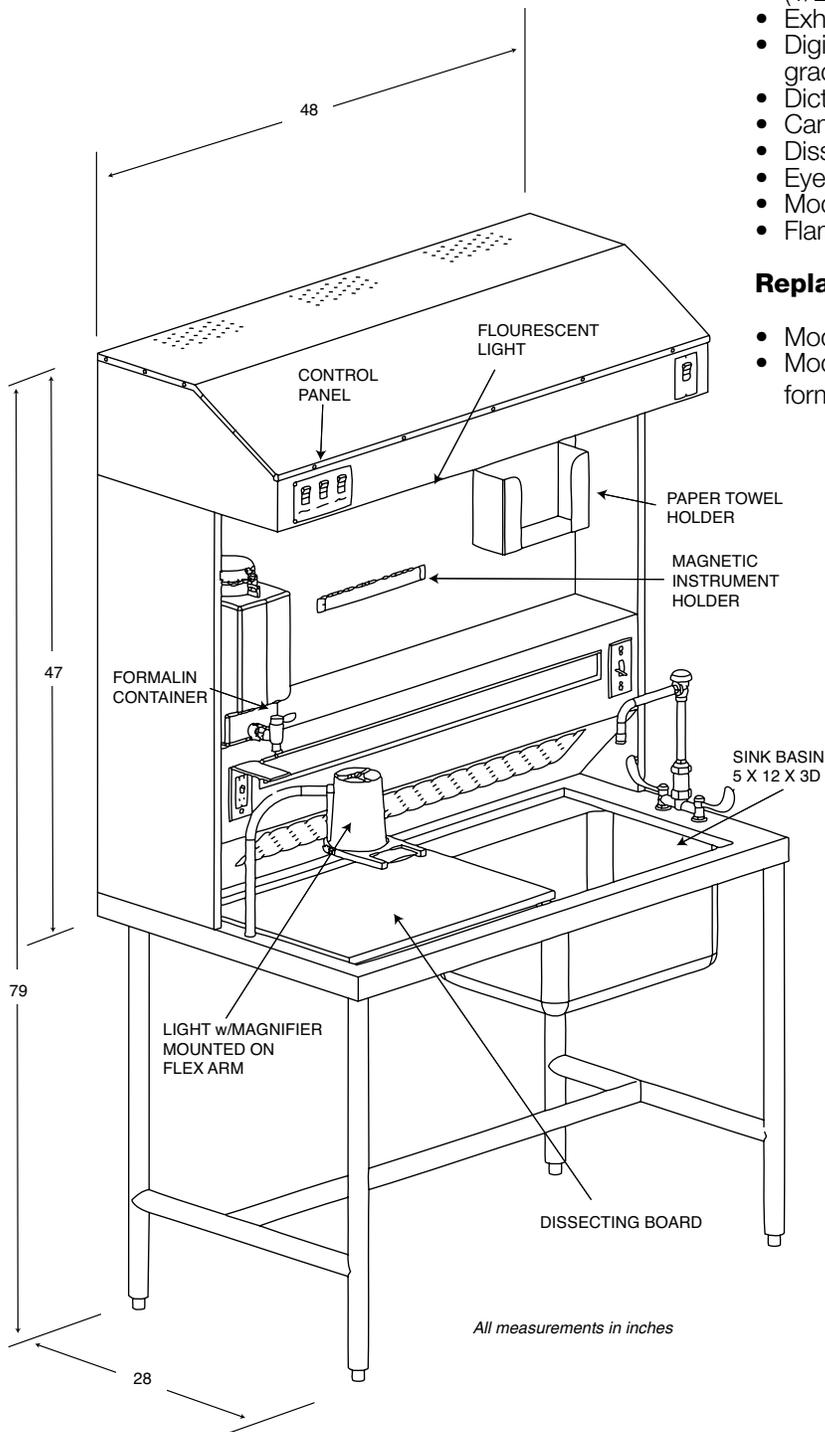

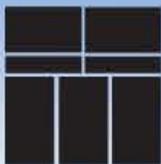

### GROSSING STATION ACCESSORIES



#### Optional Accessories

**Monitor and Keyboard Stand** designed for your flat panel monitors. Articulating arm allows placement in a variety of positions.

**Hands Free Controls** are available for easy On/Off control of the Table Rinse, Mixing Faucet, Waste Disposal and Lighting.

**Cassette Holders** conveniently holds cassettes within arms reach. Removable boxes mount to rail on front of unit.

**Flex Arm Halogen or Incandescent Lights** mounted on station to direct focused light over work area. Incandescent lights available with magnifiers. Up to two units can be placed on each station. Reticulated magnifier with fluorescent lighting allows for positioning in various areas.




**Stainless Steel C-Fold Paper Towel Holder** is mounted above the sink area


**Plexiglass Splash Shield** can be mounted over the work area to protect the user while still providing excellent visibility of the work area. Flex arm allows easy adjustability.

**Form Holder** can be mounted for storage of documents up away from any water and risk of damage.

**Three Side Rinse or End Rinse** provides a constant flow of water to the work surface to wash debris toward the sink basin. Unit can be manual valve or hands free controlled.

**Ventilated Trash Container** made of stainless steel can be provided for waste disposal. Top lid removes to hold a trash liner. 1.5" duct attaches to station exhaust system for fume removal.

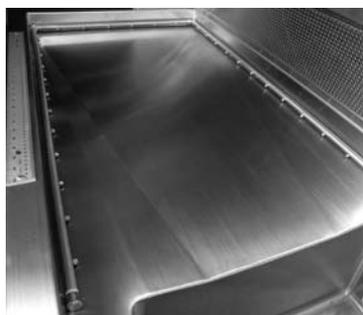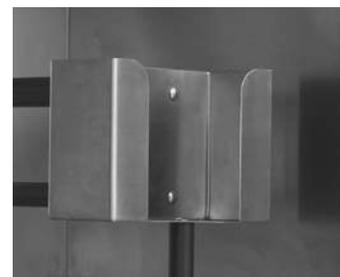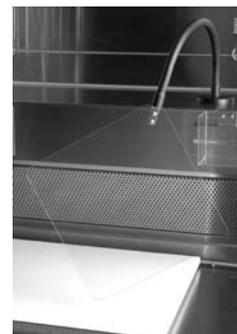

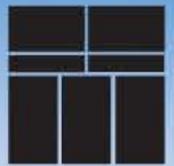

## GROSSING STATION ACCESSORIES



### Optional Accessories

**Heavy Duty** commercial disposal (1/2 hp, 3/4 hp, 1 hp)

**Sliding Access Door** for easy filter access. Filters can be changed without tools from front of station.

**Replacement Filters Model GL112** Activated Charcoal Filter for odorous fumes. **Model GL113** Potassium Permanganate Filter for formaldehyde fumes.

**Exhaust duct** connection to building ventilation system

**Base Leg Assembly GL120** 30" deep, built to customer specifications for width and length. Adjustable leveling feet.

**Base Assembly with Cabinet Model GL125** 30" deep, built to customer specifications for width and length. Adjustable leveling feet.

**Formalin Dispensing and Collecting System** Formalin disposal cup sink connected to carboy collection bottle and formalin dispensing faucet connected to carboy with dispensing pump. 2.5 gallon formalin container with electronic eye to signal alarm when supply is low. Swing door allows access to containers for replacement.

**Camera Mount** allows for use of a digital or 35mm camera.

**Dissecting Area Rinse Assembly**

**Eyewash Assembly**

**Removable Measuring Rule** in Cm/In

**Economy Grade or Stainless Steel** Digital Platform Scale (2 Kg).

**Other Optional Features Include:** Foot Operated or Voice Activated Dictation Equipment, Microphone on Flex Arm, Flammable Storage Cabinet, Video Camera Arm, Pull Out Writing Platform, X-ray Illuminator and Utility Drawer Assembly

