



MORTECH MODEL 0224-01 AUTOPSY SAW



Autopsy Saw

State of the art high-speed oscillating autopsy saw with features and technology that have been long overdue in the autopsy/pathology profession. The Mortech Autopsy Saw is an essential instrument for pathologists. The high speed oscillating action quickly cuts through bone without harming soft tissue. The saw may be used for removing the skull cap, making linear cuts, or taking small bone specimens.

Outstanding Features:

- Completely enclosed motor – use cleaning solution on damp cloth for quick and easy sterilization
- Light weight – 2 pound 6 ounces
- Easy to handle – ergonomic shaped to fit the hand for a more natural and comfortable fit
- Heavy-duty waterproof cord and plug
- Custom designed motor for added power
- Convenient watertight on/off switch
- Optional adapter available for use with M-Pact Bone Dust Collector Model 0224-06

Unit ships complete with 2.5 inch stainless steel round, universal arbor and arbor wrench



MODEL 600025 HYDRAULIC AUTOPSY CARRIER



Autopsy Saw

Model No. 0295-270

- Cool operation gives greater user comfort
- Reduced noise level
- Teflon housing cleans easily
- Ten foot cord for greater mobility
- 11 inches long, 3 lbs 4 oz

Unit ships complete with 2.5 inch diameter sectioned blade with arbor, wrench, owner's manual and registration card.

Model No. 0295-271

230 volt Euro style plug

Bone Dust Vacuum


Model No. 0295-420

- 115 volt, 8.5 amps
- Reduced noise level
- Powerful enough for even the smallest particles
- Disposable filter cartridge traps more dust particles without clogging and is easy and safe to replace
- 10 foot power cord
- 23 inches long x 15 inches wide x 23.5 inches high

Unit ships complete with Vacuum nozzle adapter and owner's manual.

Model No. 0295-410

Disposable Hepa Filter Cartridge

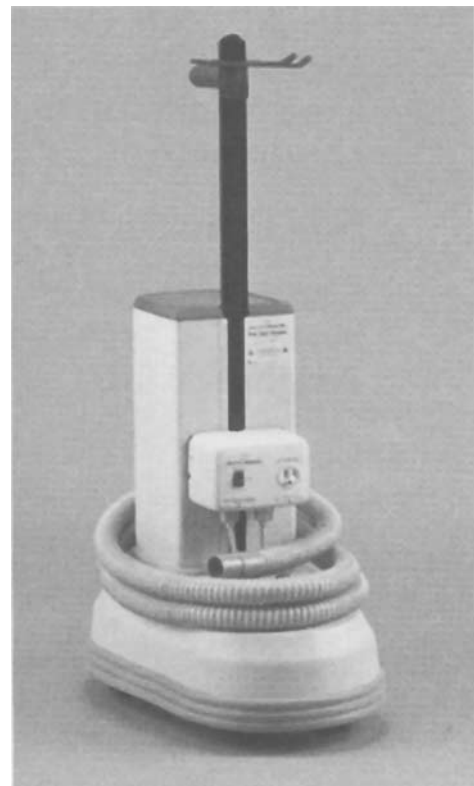



AUTOPSY SAW BLADES



BD0224-02



BDA0224-50



BD0224-03



BDA0224-51



BD0224-52



BDA0224-54



SA0224-04

Autopsy Saw Blades

Description

- 2.5 inch Round Blade without arbor
- 2.5 inch Section Blade without arbor
- Standard Saw Arbor
- 2.5 inch Round Blade with arbor attached
- 2.5 inch Section Blade with arbor attached
- 3.875 inch Spinal Column Saw without arbor
- 3.875 inch Spinal Column Saw with arbor attached

Order No

- BD0224-02
- BD0224-03
- SA0224-04
- BDA0224-50
- BDA0224-51
- BD0224-52
- BDA0224-54