

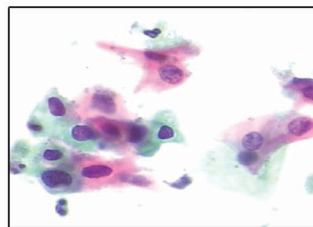


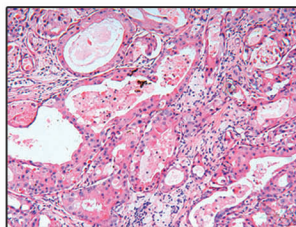

LST94 Tissue Stainer

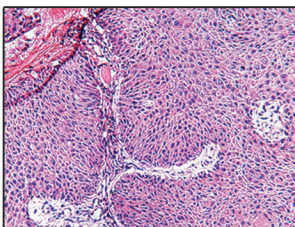

LST94 Tissue Stainer

Main characteristics:

- Protection from getting stuck: if the hanging basket comes into contact with any obstacles when entering the liquid tank during operation, it will return to the first tank so as to protect the tissue from being discarded.
- Self-adapted cover used with the multiple gasket seal produces excellent sealing performance, which minimizes the volatilization loss of the liquid.
- The organic glass housing, sealing strips and double layer active carbon filter ensure that the air in the room is not polluted.
- The larger screen liquid crystal display enables the operator to see clearly the operation condition.
- Manual operation can be used at will during running to adjust the staining effect.
- Tissue sections are washed with circulating water to ensure clean sections without contamination.

Technical Parameters:

- Quantity of liquid tank: 11
- Dimension of liquid tank: 1300ml
- Quantity of circulating water tank: 1
- Loading capacity at a time: 94 pieces of sheet glass
- Retention time in each tank: 1 second~99 minutes and 59 seconds
- Moving up and down: once every 15 seconds
- Power supply: 220V \pm 10% at 50/60Hz or 110V \pm 10% at 50/60Hz
- Power: 120W
- Outside dimensions(mm): 1200 \times 530 \times 590
- Weight: 98kg